

ABSORBABLE SUTURES



P.G.A.

Braided, violet, synthetic absorbable polyglycolic acid surgical suture, PGA-PCI-copolymer coated.

Polyglykosaure chir. Nahtmaterial. Violett, geflochten, PGA-PCI beschichtet, resorbierbar.

Tensile Strength:
65% in the 2nd week
35% in the 3rd week

Absorption:
60 - 90 days




NEOSORB (PGLA)

Braided, violet, synthetic absorbable polyglactin 910 surgical suture, coated.

Polyglactin 910 chir. Violett, geflochten Nahtmaterial resorbierbar, beschichtet.

Tensile Strength:
50% in the 3rd week

Absorption:
56 - 70 days

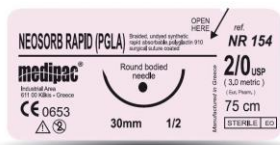


NEOSORB RAPID (PGLA)

Braided, undyed synthetic rapid absorbable, polyglactin 910 surgical suture, coated.

Polyglactin 910 chir. Rapid ungefärbt geflochten, Nahtmaterial, resorbierbar, beschichtet.

Tensile Strength:
50% in the 5th day
Total lost in 10 - 14 days

Absorption:
42 days

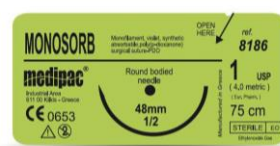


MONOSORB

Monofilament, violet, synthetic absorbable, poly (p - dioxanone) surgical suture – POO Poly-p-dioxanone chir. Nahtmaterial. Violett, monofil, resorbierbar.

Tensile Strength:
65% - 70% in the 4th week
50% - 60% in the 5th week

Absorption:
180 - 220 days

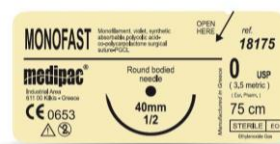

MONOFAST

Monofilament, violet, synthetic absorbable, polyglycolic acid-co-polycaprolactone surgical suture- PGCL.


Polyglykosaure-co-polycaprolactone chir Nahtmaterial. Violett, monofil, resorbierbar.

Tensile Strength:
68% - 78% in the 1st week
30% in the 2nd week
Total lost in 28 days

Absorption:
90 - 110 days

NON – ABSORBABLE SUTURES


PROPYLEN

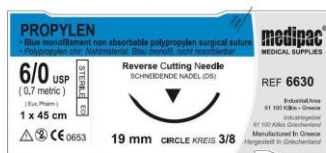

Blue monofilament non absorbable polypropylen surgical suture.

Polypropylen chir. Nahtmaterial. Blau monofil, nicht resorbierbar.


POLYAMID

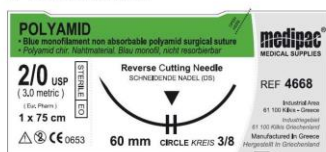

Blue monofilament non absorbable polyamid surgical suture (nylon).

Polyamid chir. Nahtmaterial.
Blau, monofil, nicht resorbierbar.

POLYESTER

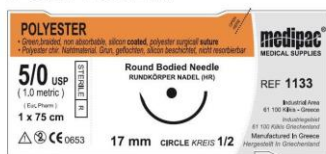

Green, braided, non absorbable, silicone coated, polyester surgical suture

Polyester chir. Nahtmaterial.

Grün, geflochten, Silikon beschichtet, nicht resorbierbar.


SILK



Black, braided, non absorbable, wax coated, silk surgical suture.

Seide chir. Nahtmaterial. Schwarz, geflochten, Wachs beschichtet, nicht resorbierbar.

PROFIMED

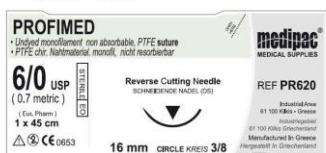

Undyed monofilament non absorbable, ptfе suture

PTFE chir. Nahtmaterial monofil, nicht resorbierbar.

SUPRAMID



Black, Pseudo-monofilament, non absorbable polyamide (nylon) sterile surgical suture

Schwarzes Nahtmaterial, monofil, nicht resorbierbar (Nylon).

For further information about details and references, please ask for our complete Medipac catalog

REY MEDICAL SA: YOUR PARTNER IN SURGERY

• **QUALITY** • **SECURITY** • **ECONOMY**

Visit our website: www.reymedical.ch or call us: 022 743 24 10