


INTRAMESH T1

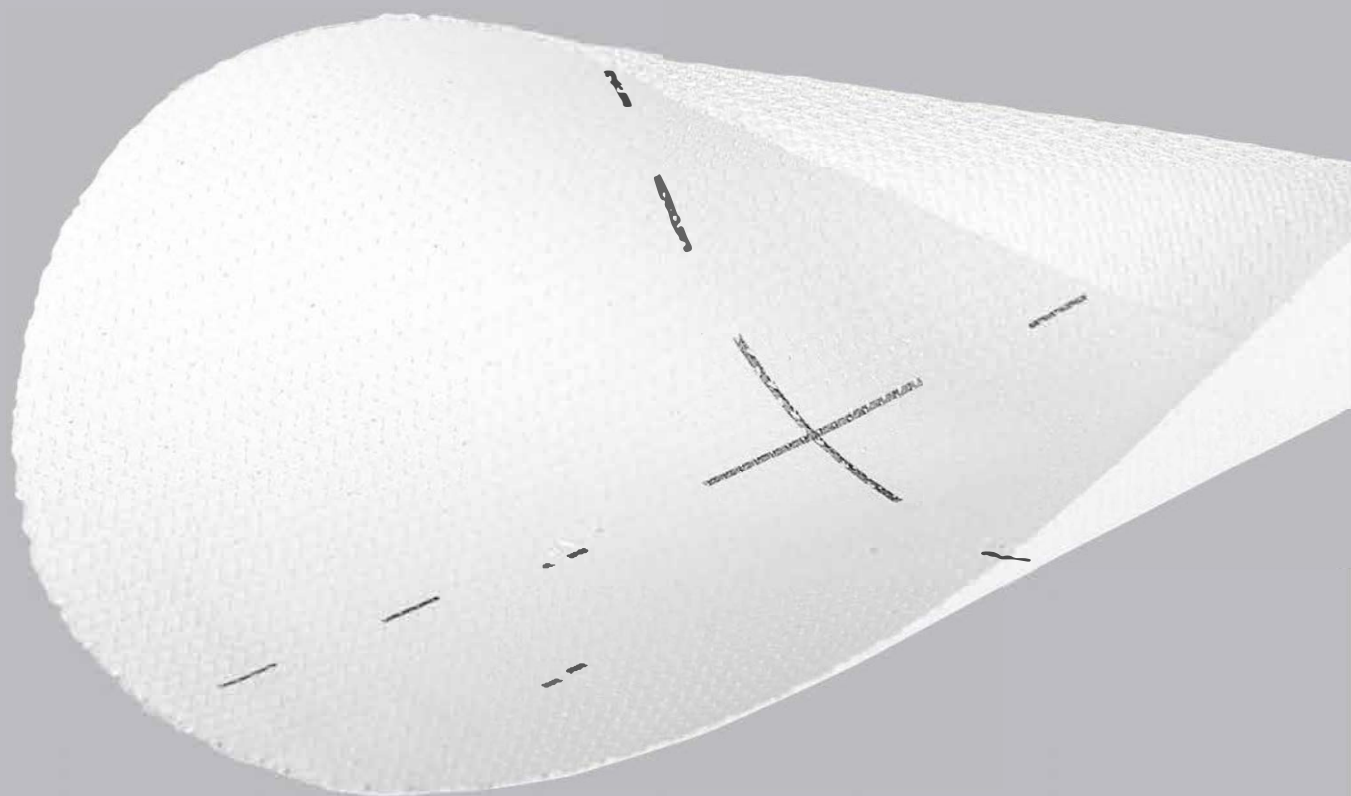

**REPAIR OF UMBILICAL
& VENTRAL HERNIAS**

**PARIETAL REINFORCEMENT
IMPLANT**

- › Safe & reliable mesh
- › Proven efficiency
- › Easy to place



INDICATIONS


Repair and parietal reinforcement for umbilical and ventral hernias

PRODUCT DESCRIPTION

These products can be used in laparoscopic surgery or open surgery. INTRAMESH® T1 products have a non-adhesive side, which is designed to avoid or reduce visceral adhesions during intra peritoneal use or extra peritoneal use in the presence of a peritoneal defect. Marking in black ink, composed of a cross and dots is on the non-adhesive side. The meshes, biocompatible and non resorbable have the advantage of having a shape memory and a great flexibility.

MATERIAL AND WEIGHT

> 120 ± 10 g/m²

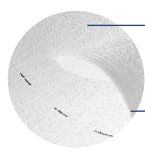

SIDE A

Smooth microperforated ePTFE
to minimize adhesions

SIDE B


Macroporous polypropylene matrix
for fast tissue ingrowth

> Proven efficiency

> Best incorporation

> Easy to handle

- Open or laparoscopic
- Thin and flexible profile
- Easy-deployment

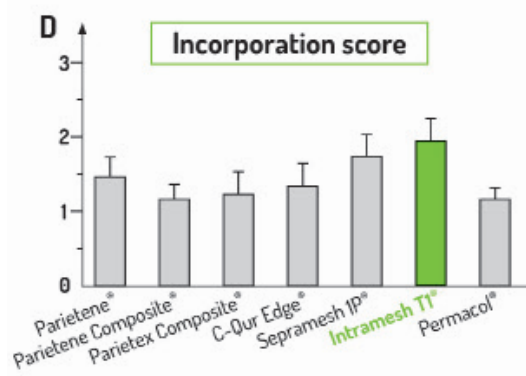

REFERENCES AND SIZES

UMBILICAL & VENTRAL

Forme	Référence	Description	Dimension (cm)
	FBIO T1 1515*	Mesh	15 x 15
	FBIO T1 2020	Mesh	20 x 20
	FBIO T1 3030*	Mesh	30 x 30
	FBIO T1 3050	Mesh	30 x 50
	FBIO T1 1015*	Mesh	10 x 15
	FBIO T1 1520*	Mesh	15 x 20
	FBBIO T1 2025*	Mesh	20 x 25
	FBIO T1 R12X*	Mesh	Ø 12

*References available in our stock at REY Médical SA

(1) Coated meshes for hernia repair provide comparable intraperitoneal adhesion prevention. Marc H.F. Schreinemacher et al. ; Surgical endoscopy (2013) 27 : 4202 -4209. Animal study.

Intramesh T1 is a class IIb medical device manufactured by COUSIN BIOTECH S.A.S. The CE conformity has been carried out by the notified body SGS Belgium NV (CE1639). The management system of COUSIN BIOTECH S.A.S is certified for compliance with ISO 13485 standard. Please read carefully the instructions for use before using the device. The IFU is available electronically at: <https://www.cousin-biotech.com/en/implant-notice>. Reference : FPMVT1GB06- 13/07/22. Non contractual pictures and texts. Specifications likely to be modified without notice. Cousin Biotech S.A.S, capital of : 340 656 € • 398 460 261 RCS Lille • N°VAT FR 34 398 460 261. Cousin Biotech is the legal manufacturer of the medical devices proposed by Cousin Surgery.

YOUR PARTNER IN SURGERY SINCE 2003

INNOVATION AND QUALITY AT THE SERVICE OF SURGERY

+41 (0)22 743 24 10 • contact@reymedical.ch • www.reymedical.ch

The best technology available on the market, at the best price, at the right time